

imp De Mayo Knee Positioner®



"Like having another pair of hands in the O.R."

***Single Locking Lever Allows for Precise Control
of Flexion, Extension, Tilt and Rotation***

impTHE OPERATIVE WORD
IN PATIENT POSITIONING


De Mayo Knee Positioner®



The De Mayo Knee Positioner® reduces surgeon fatigue and assists in the stabilization of a patient's leg during surgery. The unique *single locking lever* allows for precise control and its ease of use saves surgical time. The same lever releases the boot for testing of ligament stability and ROM testing. The positioner is easily set up in the sterile field on top of the sterile drapes.



The *single locking lever* controls flexion, extension and boot release; as well as finite internal/external rotation allowing superior access to the lateral and medial compartments.

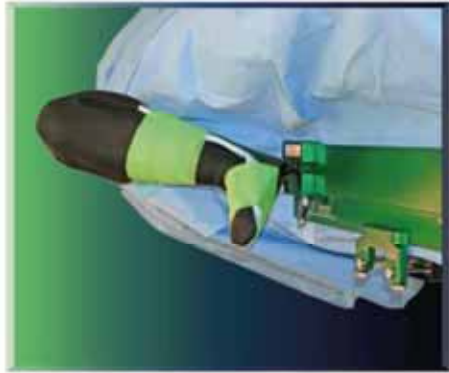

External Rotation



Internal Rotation



Single Lever Clamp securely attaches positioner to O.R. table siderail over the sterile drape without perforation. (Cat. #903)



De Mayo Push Button™ Clamp featuring Push Button release mechanism for *Ease of Use*. Significantly reduces OR set up time. (Cat. #713-302)



Flexion




Extension

imp De Mayo Knee Positioner® Accessories



5 Section OR Table Pads
(Cat. No. 398-LG Kit)

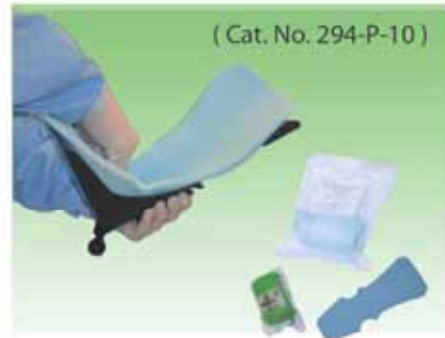

(Cat. No. 294-P-10)



Single use, sterile foam pad and cohesive wrap protect the patient from pressure ulcers, abrasions and possible neurological impairment while securing the foot into the boot.

De Mayo Universal Distractor®




(Cat. No. 907)



(Cat. No. 907-P)

Sterilization Cases Kits

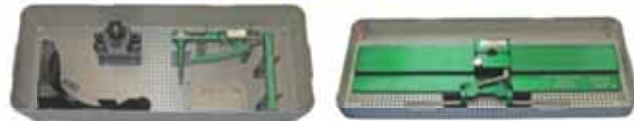

Cat. No. 706 Kit



Cat. No. 706-XS Kit

Ordering Information:

Description

- De Mayo Knee Positioner® 30" Base Plate, Aluminum Boot & Single Lever Clamp
- De Mayo Knee Positioner® 30" Base Plate, Aluminum Boot & Push Button Clamp™
- De Mayo Knee Positioner® 25" Base Plate, Aluminum Boot & Single Lever Clamp
- De Mayo Knee Positioner® 25" Base Plate, Aluminum Boot & Push Button Clamp™

Cat. No.

- 803 Kit
- 803-A713-3
- 803-S Kit
- 803-S-A713-3

Accessories:

- De Mayo Knee Positioner® Aluminum Boot with Distractor Block
- De Mayo Knee Positioner® Composite Boot with Distractor Block
- De Mayo Universal Distractor®
- Single Lever Clamp
- De Mayo Push Button Clamp™
- De Mayo Push Button Clamp™ -Long
- 5 Piece Set O.R. Table Pads for use with De Mayo Knee Positioner® (includes 3 piece O.R. Table Pads & 2 piece Armboard Pads)
- Sterilization Cases 706 Kit (includes 706-S & 706-L)
- Sterilization Cases 706-XS Kit (includes 706-S & 706-XS)

- 803-ABD
- 803-CBD
- 907
- 903
- 713-302
- 713-803
- 398-LG Kit
- 706 Kit
- 706-XS Kit

Patient Protective Pads

- De Mayo Knee Positioner® Sterile Foam Pads & Cohesive Wrap 10/Case
- De Mayo Universal Distractor® Sterile Pressure Protective Pad Set 10/Case (includes Distractor Pad, Foam Pad & Cohesive Wrap)

- 294-P-10
- 907-P



De Mayo Knee Positioner® is a registered trademark of Innovative Medical Products, Inc.
De Mayo Knee Positioner® U.S. Patent # 7,390,299
De Mayo Universal Distractor® is a registered trademark of Innovative Medical Products, Inc.
De Mayo Universal Distractor™ U.S. Patent # 8,048,082
Single Lever Clamp- U.S. Patent No. 7,003,827
De Mayo Push Button Clamp™ is a trademark of Innovative Medical Products, Inc. -Patent Pending For
© January 2015 Innovative Medical Products, Inc. All Rights Reserved



imp Innovative
Medical Products, Inc.

87 Spring Lane, PO Box 8028
Plainville, CT 06062

PH: (860) 793-0391 FAX: (866) 459-1805
Toll Free PH: 1-800-467-4944

www.impmedical.com